
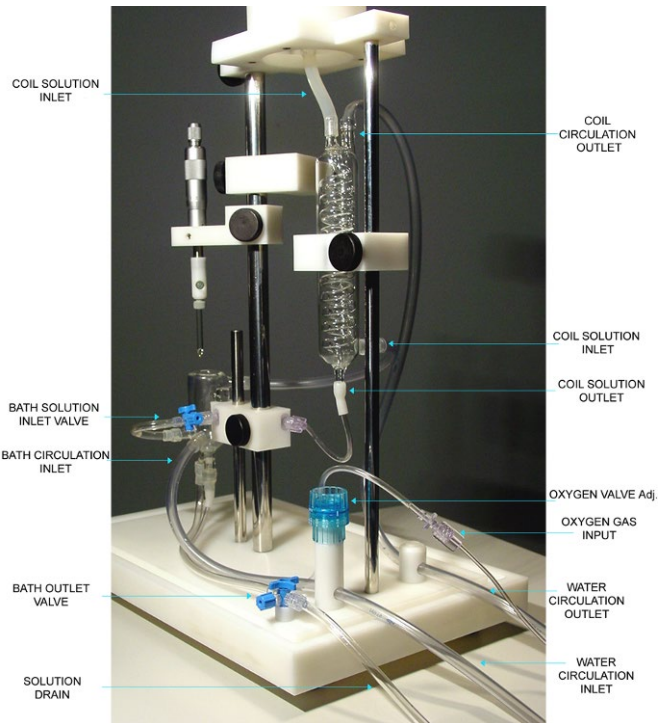


## TISSUE BATH 1, 2, 4, 8 TISSUE BATH STATIONS

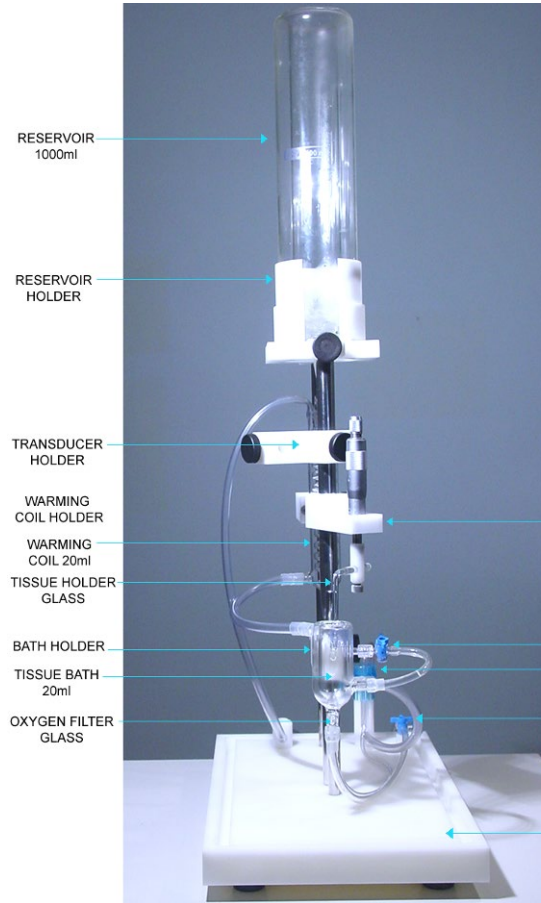

The Tissue Bath Station is completely modular and can be purchased in multiples of one unit. The System includes all the glassware, tubing, reservoir, tissue hooks and mounting accessories, force transducer and micrometer tension adjuster.

The ergonomic design of the station allows the tissue bath to be lowered away from the tissue holder so that mounting of the tissue preparation is very easy. The taps for filling and draining the bath are mounted on the tubing to avoid the risk of accidental bath breakage. The entire station is mounted on a convenient base stand, which creates a sturdy platform for the experiment. The unique design makes it easy to add or remove stations to provide the optimal solution for the requirements. Must specify tissue bath and heating coil size when ordering.

### TISSUE BATH SYSTEM COMPONENTS

- |   |   |
|---|---|
| 1 x Tissue Holder—stainless steel; reorder as RXHOLDER-TR             | 1 x Circulator pump—15 W; 500 ml/min  |
| 1 x Tissue Holder—glass; reorder as RXHOLDER-G                        | 1 x Warming Coil; reorder as RXCOIL   |
| 1 x Tissue Holder —stainless steel; reorder as RXHOLDER-S             | 1 x Oxygen Filter; reorder as RX02FILTER  |
| 2 x Tissue Clip—stainless steel; reorder as RXCLIP                    | 1 x Bath —reorder as RXBATHsize for 5,10, 20, or 30 ml  |
| 2 x Triangle Tissue Clip—stainless steel; reorder as RXCLIP-TRI       | 1 x Reservoir—1000 ml; reorder as RXRESERVOIR   |
| 1 x Replacement Warming Reservoir 400 ml; reorder as RXWARMING        | Mount Accessories Kit; reorder as RXMOUNT   |
| 1 x Integrated heater—1,600 ml volume, programmable temp. 20° - 44° C | Field Stimulation Electrode; reorder as STIMHOLDER for STM100C, BSLSTIMHLD for BSL Stimulator |
|   | 1 x Micrometer-transducer assembly  |
|   | 1 x 3-way Rotary Valve  |
|   | 1 x Power Supply – 110V/60 Hz or 220V/50 Hz   |

**See also:** Tissue Bath Accessories next page  
BIOPAC Circulators

## TISSUE BATH ACCESSORIES / REORDER PARTS

### Tissue Holders and Clips



RXHOLDER-S



RXHOLDER-G



RXHOLDER-TR




RXCLIP



RXCLIP-TRI

### Warming Coil

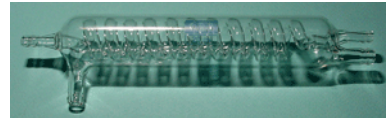

### Oxygen Filter




### Tissue Bath




### Reservoir



### Mount Accessories

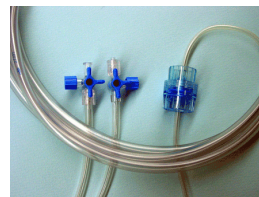

### Field Stimulation Electrode

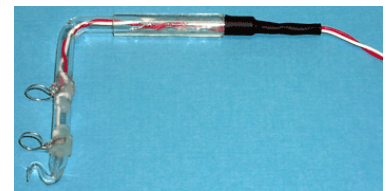

<b>RXHOLDER-S</b>	Tissue Holder (stainless steel) 77.34 mm high x 55.06 mm wide – Qty 1
<b>RXHOLDER-G</b>	Tissue Holder (glass) 67.46 mm high x 57.85 mm wide – Qty 1
<b>RXHOLDER-TR</b>	Triangle Tissue Clip Holder (stainless steel) 15 mm high x 9 mm wide – Qty 2
<b>RXCLIP</b>	Tissue Clip (stainless steel) 15 mm high x 5 mm wide – Qty 2
<b>RXCLIP-TRI</b>	Triangle Tissue Clip for Rings (stainless steel) 15 mm high x 12 mm wide – Qty 2
<b>RXWARMING</b>	Replacement Warming Reservoir 400 ml
<b>RXCOIL</b>	Warming Coil
<b>RXO2FILTER</b>	Oxygen Filter (glass)
<b>RXBATHsize</b>	Tissue Bath: RXBATH5 (5 ml,) RXBATH10 (10 ml,) RXBATH20 (20 ml), RXBATH30 (30 ml)
<b>RXRESERVOIR</b>	Reservoir 1000 ml
<b>RXMOUNT</b>	Mount Accessories Kit
<b>STIMHOLDER</b>	Field Stimulation Electrode for use with STM100C
<b>BSLSTIMHLD</b>	Field Stimulation Electrode with BNC cable termination for use with BSL Stimulator